


GladiatorFit

GladiatorFit Massage Gun


gladiatorfit.ch

Product Structure Chart

GladiatorFit

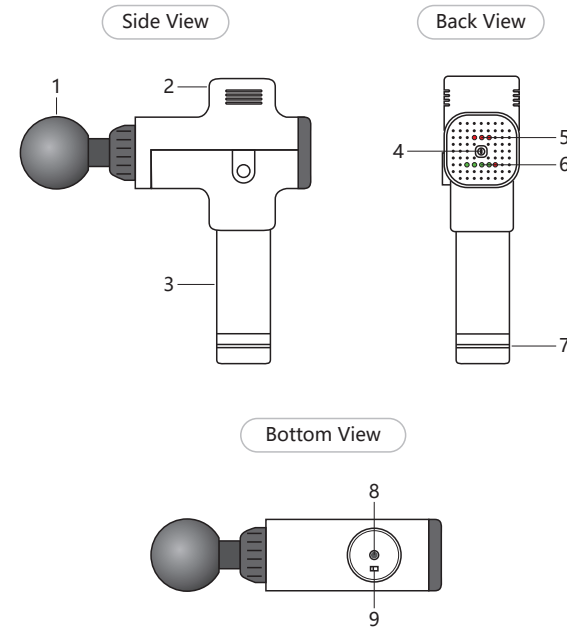

- Globule head**
Location of muscle groups
- Forked head**
Suitable for vertebral position
- Flat head**
Applicable to all parts of the body
- Spiral head**
Joint parts of the body

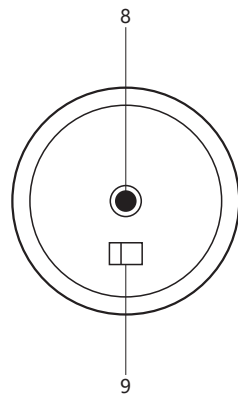

- 1, Removable, replaceable applicators. Pull firmly to remove
- 2, High speed brushless motor
- 3, ACwall charter(24VDC-1A)
- 4, Power switch
- 5, Frequency level indicator LED(blue)
- 6, Battery level indicator LED(green, red)
- 7, Battery release button
- 8, DC charging port
9. Power ON/OFF button

SPECIFICATIONS

AC Charger: 110-240V 50/60Hz-.7A-24VDC-1A

Pulse(BPM): Level 1-1800, Level 2-2400, Level 3-3200

Weight: 2.5 lbs/1.1 kg (3.0 lbs/ 1. 4 kg including charger)

Power device on by moving the power switch on the bottom of the battery to the ON position Power device off by moving the power switch to the OFF position.

With power switch in ON position, press the power button once to turn device on to level one, twice to turn device to level two and three times for level three. In the same way, you can adjust to level four or level five. Pressing the power button again will turn the device off. The blue LED lights will indicate the corresponding power level.

Please turn off the switch after the massage is completed.

Massage desired body part at the desired speed applying pressure as needed.

The product has a built-in protection device that automatically stops when it is running for 10 minutes.

Charging



- Fully charge battery for up to six hours before first use either attached or detached from the device.
- To charge, connect the DC end of the supplied 24V adapter to the battery's charging port and plug AC side cable into wall outlet.
- Battery may be recharged at any time and at any battery level.
- It is not recommended to fully deplete the battery to the red LED level.
- The average run-time is over two hours depending on speed level and pressure applied during use.
- Never leave the device unattended whilst charging, and always charge on a non-flammable surface.

SAFETY INSTRUCTIONS

DANGER:

TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS.

Warnings and Caution



1. For adult use only.
2. Use only on a dry clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region.
3. Use the massage device only on the body's soft tissue as desired without producing pain or discomfort. Do not use on head or any hard or bony area of the body.
4. Use only the applicators that provide the best desired effect.
5. Bruising can occur regardless of control setting or pressure applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.
6. Keep fingers, hair or any other body part away from the shaft and backside of the applicator as pinching may occur.
7. Do not place any objects into the ventilation holes of the massage device.
8. Do not immerse in water or allow water to enter the ventilation holes of the massage device.
9. Do not drop or misuse the massage device.
10. Only recharge with supplied charger.
11. Carefully examine device and battery before each use.
12. Do not tamper or alter massage device in any way.
13. Never leave the massage device operating or charging unattended.

Maintenance



Use a slightly damp towel to clean device and a dry with a soft cloth

For storage or travel, remove battery and store in supplied case.

Warnings and Caution

Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.

Do not remove screws or attempt to disassemble.

Do not operate continuously for more than one hour. Allow device to rest for 30 minutes before re-using.

Unplug the charger from device battery after charging or prior to use.

For adult use only. Do not use if injured. Consult your doctor before using this product.

Use product only according to Instructional Use materials available at Product instructions attachments.

Warnings and Caution



PLEASE DO NOT USE THE MASSAGE DEVICE, OR ANY VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY:

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device but we advise you to consult a doctor first.

Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above. Practical experience has shown that integration of vibration exercise into a treatment plan is advisory in a number of cases. This must be done on the advice of and in accompaniment with a doctor, specialist or physiotherapist.

Warranty Statement



The equipment is manufactured from the time of purchase. If the product fails due to defects in materials or design process within one year, the company will repair or replace parts or replace new products for free, except in the following Cases:

1. Improper use of human or damage to the equipment caused by transportation.
2. Unauthorized disassembly and repair of this equipment.
3. Failure to follow the instructions.
4. The product is damaged due to the abnormal storage or maintenance environment of the customer.
5. If the proof of purchase date is not provided, the company will reserve the right to refuse the warranty.